

Biopsy?


the complete solution

Because of the nature of the sample, the laboratory has to take great care that a biopsy is well protected during tissue processing.

This is where CellPath USA can help the histotechnologist. We offer a range of biopsy handling products that allow you to choose the right option for your tissue.

- CellSafe
- CellSafe5
- Biopsy Sponges
- Biopsy Bags
- Tissuewrap

Whatever your need - we have the solution!

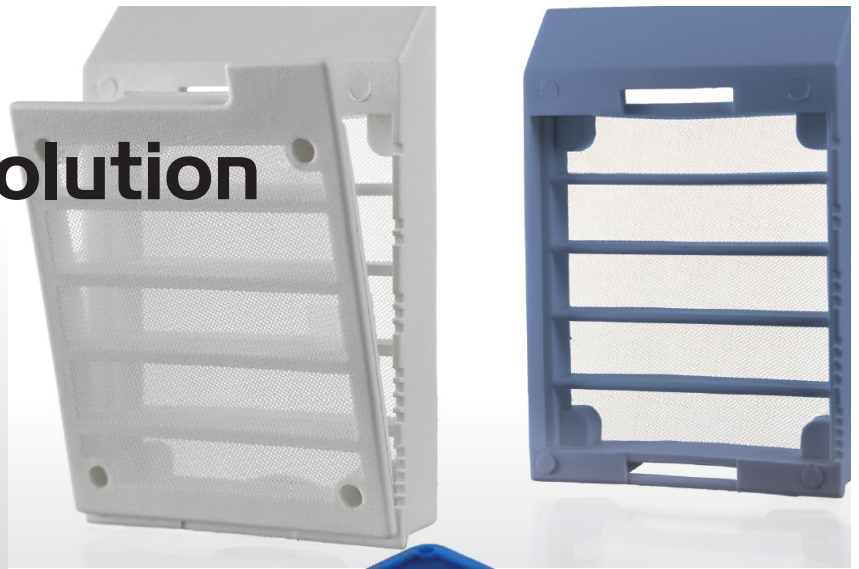

BioPad Biopsy Sponge

CellPath's BioPad biopsy sponges are cut to a size that fits all standard processing cassettes. BioPads allow fluid to transfer thoroughly and safely around the tissue sample while holding the tissue specimen safely in place.

- Cut to size to fit most major processing cassettes
- Thin smooth section to minimize imprinting on specimens
- Prevents loss of small specimens
- Low cost
- Efficient
- Available in black or blue sponge

Part No	Description	Qty
EBA-0101-03A	BIOPAD BLACK	1000
EBA-0106-03A	BIOPAD BLUE	1000

CellSafe+ Biopsy specimen handling

The CellSafe+ consists of two hinged interlocking frames with an integral nylon mesh. This extra fine 200µm mesh ensures a high level of specimen security, zero tissue damage and low reagent carry-over.

- New Blue Mesh provides better contrast between the biopsy and the CellSafe
- X-Ray transparent mesh for breast work
- Extra fine mesh provides up to 20 times less carry-over of solvents than sponges
- Zero tissue damage or artifacts associated with sponges
- Substantial time saving compared to folding biopsy in tissue paper
- Low reagent carry-over.

Part No	Description	Qty
EBE-0201-02A	CELLSAFE+ BIOPSY CAPSULE WHITE	100
EBE-0201-03A	CELLSAFE+ BIOPSY CAPSULE WHITE	1000
EBE-0301-02A	CELLSAFE+ BIOPSY CAPSULE BLUE	100
EBE-0301-03A	CELLSAFE+ BIOPSY CAPSULE BLUE	1000

CellSafe 5

With increasing pressure to process larger numbers of small biopsy specimens, CellSafe 5 reduces the pressure on laboratory staff by enabling them to perform five core biopsies in numbered wells, which can be sectioned together for rapid and cost effective preparation.

CellSafe 5 is a new five well cassette developed for the simultaneous preparation of core biopsies, specifically breast and prostate – two of the most costly diagnostic samples to deal with. Each well is numbered for identification and, when combined with a robust protocol will allow fast, accurate and cost effective preparation.

The cassette is available in two colors (white and blue) and features an integral X-ray lucent fine mesh for more efficient processing. The blue cassette has been produced to aid contrast and visualization of white or light colored biopsies, normally experienced with breast tissue.

Part No	Description	Qty
EBG-0201-02A	CELLSAFE 5 BIOPSY CASSETTE WHITE	100
EBG-0206-02A	CELLSAFE 5 BIOPSY CASSETTE BLUE	100



System III Hex Biopsy Cassette

The System III Hex Biopsy Cassette differs from conventional round hole biopsy cassettes in that it has hexagonal instead of round perforations arranged in a honeycombed structure to better utilize the available surface area. This unique arrangement increases the open area by approximately 13 percent thereby increasing reagent flow-through and improving specimen staining and analysis time.

A potential cost saving benefit is that there is less carry-over of chemicals, so reducing the frequency of chemical change.

Part No	Description	Qty
EAJ-0102-10A	BIOPSY CASSETTE WHITE	500

The System III Biopsy Cassette is also available in red, blue, grey, green, orange, turquoise yellow, peach and brown.

TissueWrap

TissueWrap Biopsy Paper is specifically designed for processing biopsy specimens. The tissue paper is cut to a convenient ready-to-use size. This non-absorbent paper allows fluids to flow through freely while retaining small tissue fragments. TissueWrap is supplied in a dispenser pack for increased convenience.

Part No	Description	Qty
EBB-0102-10A	TISSUEWRAP	500



Call Toll Free on **1-855-CELPATH**
Email us at **sales@cellpath.us.com**
Buy Online at **cellpath.us.com**

CellPath 
INNOVATION IN CELLULAR PATHOLOGY