

# Sample Protection at Any Temperature




0 °F      20 °F      40 °F      60 °F      80 °F      100 °F      120 °F

Inadequate storage conditions commonly leave patient samples unusable; leading to costly redraws, frustrated physicians, and unhappy patients. The SampleSafe is designed to ensure that sensitive specimens reach the laboratory as fresh as when they were collected. By selecting materials that minimize energy transfer from the surroundings, and by packing over 1" of insulation into the walls, base, and lid, the SampleSafe is designed to protect samples from the most extreme environmental conditions.

## FEATURES

1. Unmatched air or foam insulation properties guarantee the highest sample integrity possible
2. UV and moisture resistant molded design provides years of quality service
3. Separate compartment for cold packs prevent unwanted contact with patient samples
4. Folding carry handle simplifies transport inside or outside
5. Recessed lock protected from ice build up
6. Optional mounting bracket for securing the SampleSafe
7. False bottom isolates cold packs from sample

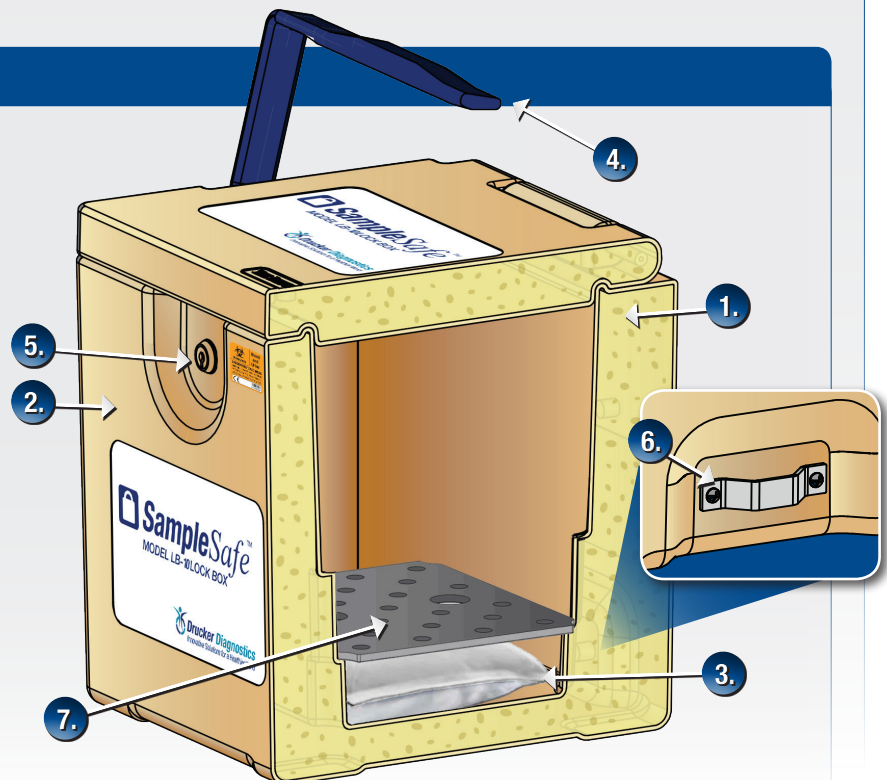


(Shown with foam insulation and key lock)

## PERFORMANCE

The superior thermal characteristics of the SampleSafe excel in maintaining sample storage temperatures in a variety of different ambient conditions. When used in combination with the optional cooling packs, sample temperatures within the SampleSafe can be safely maintained for up to 10 hours.

With the ability to maintain storage temperatures nearly 14 times longer than an insulated steel lock box, the SampleSafe ensures that maximum sample integrity is preserved (See the graph to the right).

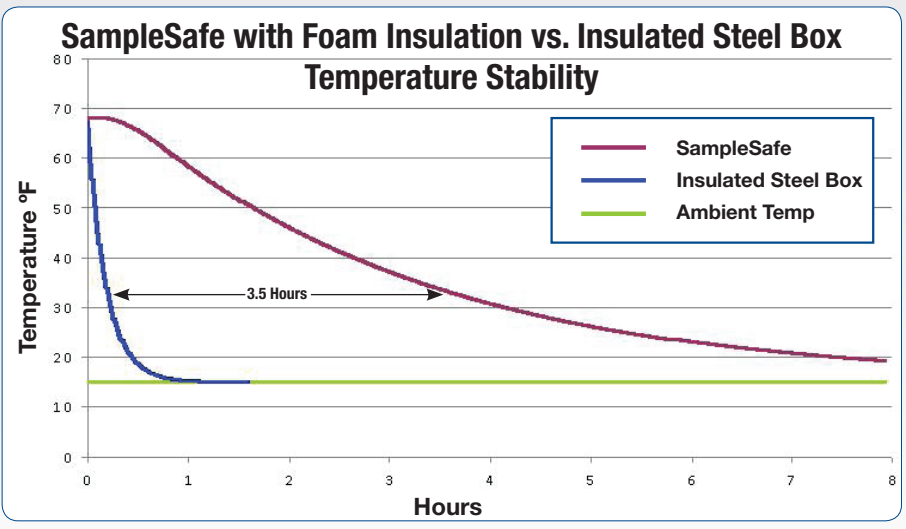

## ORDERING INFORMATION



**SampleSafe Model LB-1**  
 The SampleSafe Model LB-1 is air insulated to provide superior sample protection. Included with the Model LB-1 is a convenient carrying handle and a keyed lock. A combination lock (Catalog No. 22029579) can be substituted for the keyed lock for nominal charge. A wide variety of accessories are available to customize the LB-1.



**SampleSafe Model LB-10**  
 The SampleSafe Model LB-10 is foam insulated to provide the ultimate in sample protection. Included with the Model LB-10 is a convenient carrying handle and a keyed lock. A combination lock (Catalog No. 22029579) can be substituted for the keyed lock for nominal charge. A wide variety of accessories are available to customize the LB-10.

## SPECIFICATIONS

**SampleSafe Model LB-1**

Overall Dimensions	Height	13.5 in.
	Width	12 in.
	Depth	12.5 in.
Weight		6.7 lbs.
Internal Volume		868 in <sup>3</sup>

**SampleSafe Model LB-10**

Overall Dimensions	Height	13.5 in.
	Width	12 in.
	Depth	12.5 in.
Weight		8.5 lbs.
Internal Volume		868 in <sup>3</sup>

## OPTIONAL ACCESSORIES

### Replacement Keyed Lock

Includes 2 keys and lock to securely lock samples in the SampleSafe Lock Box.

### Replacement Key

Replacement key used with key lock (Catalog No. 22029577).

### Combination Lock

3 digit combination dial lock capable of setting 1000 different possible combinations.

### False Bottom

Rigid plastic plate used to separate samples from the cold packs (Catalog No. 22029585) preventing unwanted ice pack contact with sensitive samples.

### Under Door Bracket Kit

Includes reinforced under door bracket and hardware used to secure the SampleSafe to interior or exterior doors. Compatible with all doors up to 2 inches thick.

### Anchor Bracket Kit

Includes metal bracket, two tamper proof screws and Torx key. Allows the SampleSafe to be secured using anchor cable (Catalog No. 22029584).

### 6 Foot Anchor Cable

6 foot long x ¼ inch diameter vinyl coated metal twined anchor cable used with the Anchor Bracket kit (Catalog No. 22029582).

### Rubber Feet Kit

Includes 4 thermally insulated rubber feet with mounting hardware to isolate the lock box from the floor or pavement.

### Cold Pack

6 inches wide x 9 inches deep x 0.8 inches thick cold pack used to aid in preserving sample integrity in warm environments. Fully tested to provide ideal storage for up to 10 hours.

### One Color Screen Printing

A one color screen print of a company name and/or logo can be applied to the SampleSafe. Contact your Fisher Healthcare Account Manager for details.



Catalog No. 22029577



Catalog No. 22029578



Catalog No. 22029579



Catalog No. 22029580



Catalog No. 22029581




Catalog No. 22029582



Catalog No. 22029584



Catalog No. 22029583

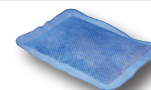

Catalog No. 22029585



Your Name Here

© 2014 Thermo Fisher Scientific Inc. All rights reserved.  
All trademarks are the property of Thermo Fisher Scientific Inc. and its subsidiaries.

#### In the United States:

For customer service, call 1-800-640-0640  
To fax an order, use 1-800-290-0290  
To order online: [www.fisherhealthcare.com](http://www.fisherhealthcare.com)

 **Fisher  
HealthCare**  
Part of Thermo Fisher Scientific